

# Piezo Cross Stage | PPX-32CR NEW

www.micronixusa.com | info@micronixusa.com | phone: +1 (949) 480 0538 | fax: +1 (949) 480 0538

The PPX-32CR is a low profile, integrated XY stage with a height of only 14 mm. Crossed roller bearings guarantee smooth and stable motion for loads up to 10 N. It utilizes a multi-phase piezo motor resulting in high speed (> 10 mm/s) and high blocking force (> 2 N). The PPX-32CR is available in open loop or with a linear encoder. Closed-loop encoder resolution of 2 nm is achievable. Versions capable of operation in vacuum ( $10^{-9}$  mbar) are available. The PPX-32CR is compatible with the MMC-100 and MMC-110 controllers.

## KEY FEATURES

- Travel range of 18 x 18 mm
- Open aperture
- 2 nm closed loop encoder resolution
- Load capacity up to 1 kg
- Crossed roller bearings
- Low profile, 14 mm height
- Vacuum versions available

## TECHNICAL DATA






<b>Travel range [mm]</b>	<b>18 x 18</b>		
<b>Straightness / Flatness [<math>\mu\text{m}</math>]</b>	$\pm 2.5$		
<b>Pitch [<math>\mu\text{rad}</math>]</b>	$\pm 200$		
<b>Yaw [<math>\mu\text{rad}</math>]</b>	$\pm 100$		
<b>Weight [g]</b>	40		
<b>Motor option</b>	<b>Piezo Motor</b>		
<b>Speed, max [mm/s]</b>	2 (MMC-100), 10 (MMC-110)		
<b>Encoder option</b>	<b>None (open loop)</b>	<b>Analog (1 V<sub>pp</sub>)</b>	<b>Digital (RS-422)</b>
<b>Resolution, typical [nm]</b>	1	10	2
<b>Repeatability, bi-directional [nm]</b>	n/a	$\pm 50$	$\pm 50$
<b>Repeatability, uni-directional [nm]</b>	n/a	50	50
<b>Materials</b>	aluminum body, steel bearing (other materials i.e. stainless steel, titanium, etc. available upon request)		

## ORDERING INFORMATION

PPX-32CR- **1** **1** **0**

<b>DRIVE</b>	Piezo Motor, PM-002 .....	1
<b>TRAVEL</b>	18 x 18 mm .....	1
<b>ENCODER</b>	None .....	0
	Analog (1 V <sub>pp</sub> ) .....	2
	Digital (RS-422) .....	3
<b>LIMIT SWITCH</b>	None .....	0
	<b>ENVIRONMENT</b>	Atmospheric .....
	High Vacuum, $10^{-6}$ mbar .....	6
	Ultra High Vacuum, $10^{-9}$ mbar .....	9

Load, max	F <sub>x</sub> [N]	F <sub>y</sub> [N]	F <sub>z</sub> [N]	M <sub>x</sub> [N-m]	M <sub>y</sub> [N-m]	M <sub>z</sub> [N-m]	k <sub>ax</sub> [ $\mu\text{rad/N-m}$ ]	k <sub>ay</sub> [ $\mu\text{rad/N-m}$ ]
PM-002	2	2	10	0.6	0.6	0.6	-	-

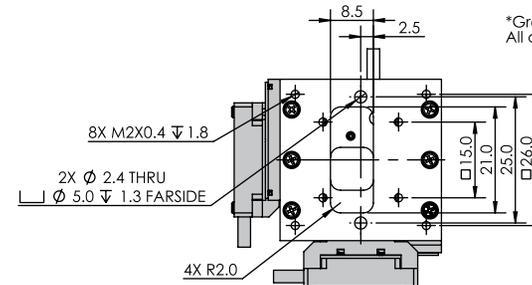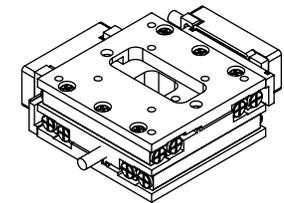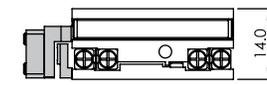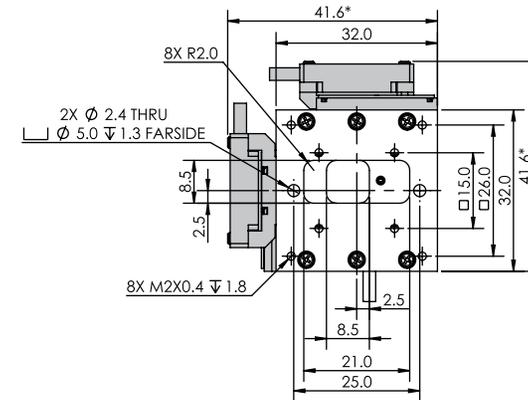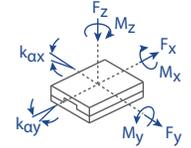

\*Grey parts are for closed loop version only  
All dimensions in millimeters

Specifications are subject to change without notice.