


SPIDER Amplified Photodetector & ADC

SPIDER is an **amplified Si-InGaAs photodetector** with **programmable gain** and an **embedded 24bit data acquisition system**. The employed two-color detector provides a straightforward beam alignment and a **wide spectral response**, ranging from **320nm to 1700nm**. The Si and InGaAs photodiodes are amplified simultaneously with two independent low-noise amplifiers, with 8 programmable gains and up to 56kHz bandwidth.

Key Features

- Amplified Detector & Analog-to-Digital Converter in a Unique Device
- Broad spectral coverage with a single-pixel element
- Fiber coupling available
- 8 programmable gains
- Up to 56 kHz bandwidth
- 2 channels, 24 bit ADC, up to 120 kSPS/channel
- 2 BNCs for analog outputs
- 3 GPIOs (general purpose input/output pins)
- USB connection to PC
- Plug&Play software and DLLs available



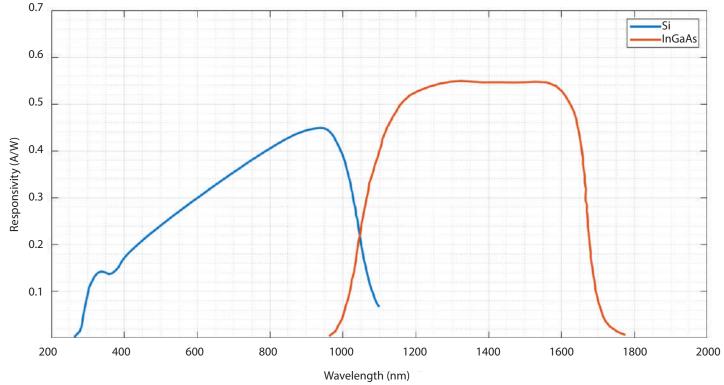



Technical specifications

Optical Specs

Detector element (Photosensitive area)	Spectral response range [µm]	Peak sensitivity wavelength λp [μm]	Responsivity R λ= λp [A/W]	I _D	Detectivity D* $\lambda = \lambda p$ [cm*Hz ^{1/2} /W]
Si (2.4 x 2.4mm)	0.32 to 1.08	0.94	0.45	130pA	1.4 x 10 ¹³
InGaAs (Ø1mm)	1.01 to 1.7	1.55	0.55	1nA	3.5 x 10 ¹²

Transimpedance Amplifier Specs

# of Programmable Gain	8							
	Selected Gain							
	#1	#2	#3	#4	#5	#6	#7	#8
Gain $[dB_{\Omega}]$	64	73	83	93	103	113	123	133
Bandwidth (-3dB) [kHz]	56	55	54	53	52	49	39	29
Output Noise* [nV/VHz]	20	22	26	36	61	121	265	680

* measured at 20 kHz

Analog to Digital Converter Specs

# of Predefined Sampling	g Rates	11 (120, 100, 80, 64, 32, 16, 10, 5, 2.5, 1, 0.5 kSPS)					
Sampling Rate [kSPS]	BW [kHz]	SNR [dB]	RMS Noise* [V]	BITs	ENOB		
120 (max.)	52	113	1.4 E-5	24	19.2		
0.5 (min.)	0.217	136.8	1.8 E-6	24	22.2		

* selected gain #=1.

for the details at all the sampling rates, please contact us.

Other Specs

	Min	Тур	Max	
Power Consumption [W]		3.5		
Power supply voltage [V]	10	12	14	
Output series resistance $[k\Omega]$	-	1.5 (on the BNC)		
Output analog voltage [V]	0	0 - 10 (on the BNC)		
Weight [g]		218		
Package size [mm]	1	101.5 x 81 x 28.5		
Operating Temp. [°C]	10		40	
Storage Temp. [°C]	-20		70	